
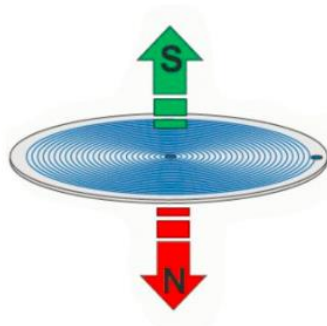
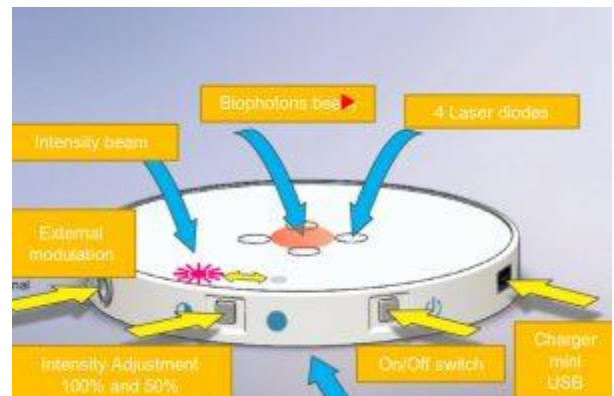


ULTRACOMPACT LASER DISC ULD 117

Amazing therapy effects using modulated soft laser and LED phototherapy, magnetic and scalar field therapy, remedies maker and water energizing all combined in the most compact therapy device.

MAIN FEATURES

- Four laser diodes, each 12 mW @785 nm for best biological effects and deep penetration; laser class 1M
- Photons at 633 nm for stimulation of collagen production
- Magnet field therapy modulated with the biological beneficial global frequency
- Tesla Scalar coil with guaranteed 40% scalar waves, 150 micro Tesla
- Modulation with Schuman Waves 7.8-8 Hz pulse for best bio-energetic information transfer and integration
- Possible modulation through external connected audio signals such as MP3 (do not use mobile phone or tables)
- New remedy imprinting, transfer of bio-signals directly to your body
- Poly-frequency modulation with solar white noise (same light spectrum as in nature)
- Miniature extra powerful and efficient rechargeable lithium battery 100-400 mA lasting for a continuous usage of 2 hours; rechargeable through USB cable
- Main unit: Dia 95 mm x 11 m height, weight 60 grams

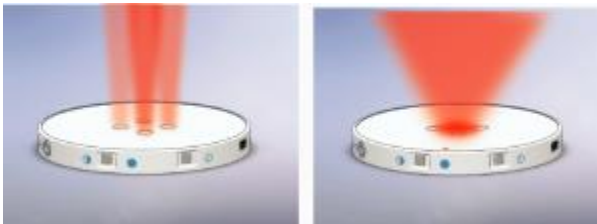

ULD 117 magnetic-tesla coil modulation

SCOPE OF DELIVERY

- Ultra compact laser disc ULD 117
- Mini USB cable
- Audio signal connecting cable
- Portable handle (ULD Mobil)
- Tesla Transmitter (USB 315)

HOW IT WORKS:

The laser diodes and LED are stimulating the water molecules inside the body (over 76%), structuring the water in hexagonal clusters (jelly water), which go into resonance with the cells. Sick cells are thus stimulated to detoxify, relieve pain and regress inflammation, boost cellular communication, repair and restore. This therapy can be applied directly on the affected body zones as well as by imprinting water and remedies which are consumed thereafter.


The magnetic field generated by a Tesla coil contains over 40% scalar waves, the frequencies carrying information; this information interacts with the DNA and microtubules, correcting and restructuring the sick body parts with healing and beneficial information. The synergy between light (photons), magnetic field together with the scalar information carriers, create a very strong healing effect. The body goes back to its coherent state of homeostasis.

ULD 117 can be used as a therapy device for all health conditions, as well as in wellness, SPA and fitness studios.

ULD 117 IS SIMPLE AND ROBUST, EASY TO OPERATE AND TO APPLY THE THERAPY


1. Switch ON/OFF to start and shutdown the ULD 117
2. Select light intensity using the keys for 100% and 50%
3. ULD 117 comes with audio port for external devices; this is a very interesting application as you can connect a MP3 device, or other devices delivering beneficial frequencies generated by music, mantras, prayers, meditation guidance, other devices
4. ULD 117 has an USB charger port for the high performant Lithium battery
5. The subtle flickering of the Laser and LED beams is the result of the modulation with the solar noise; our bodies are habituated to live in this solar modulated light, which gives the cells the "life signal"; with this light our bodies produce cholesterol sulfate in the skin, which is the building and repair block for all cells and tissues with a major function also in cellular detoxification.
6. The modulation with external audio frequencies will synchronize both the magnetic, scalar and photonic field; this increases the therapy effect by 500 times as compared with the classical PEMF, Phototherapy, scalar therapy applied separately.
7. The scalar-magnetic field is modulated with Schuman resonance (global frequency at 7.8-8 Hz); this field radiates about 2 meters from the bottom and the top of the device, penetrating your body with the south pole towards the head
8. If you remove the backside plate (unscrew the four screws), you will find 4 mini-containers; you can fill in these containers beneficial homeopathic globules, a gold nugget and close the plate thereafter; every time the device will be activated it will radiate these biosignals containing gold frequencies; the gold effect is to boost the

immune system, to stimulate your spirituality and emotional calmness, to relax but also to help better focus and concentrate. Pay attention to handle with care this operation.

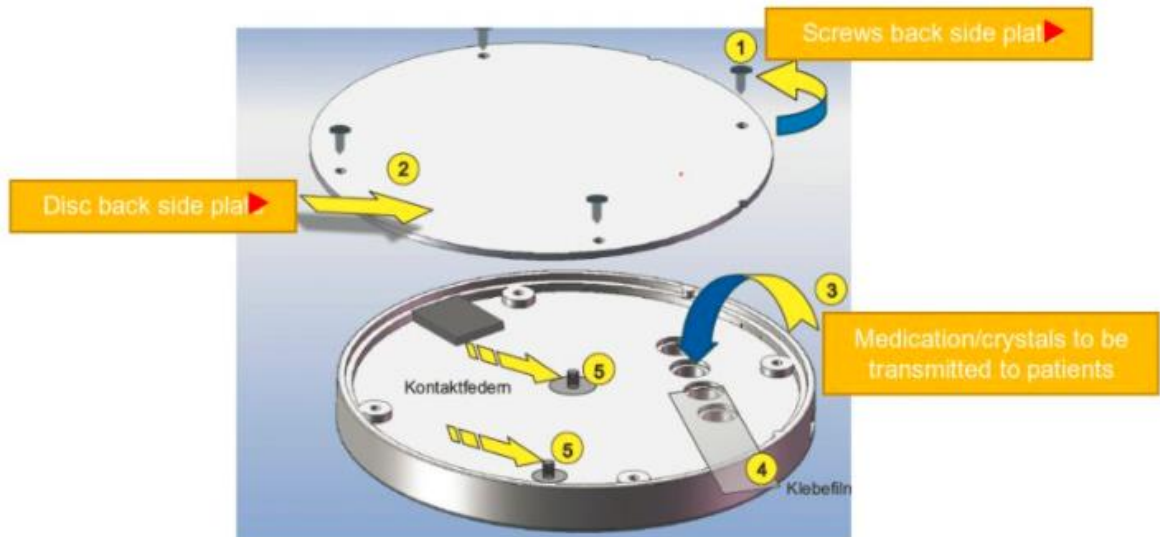

ULD 117 medication containers

9. You can imprint with ULD 117 water, other drinks, but also remedies; when you work on imprinting remedies, stay minimum 5 meters away from Wi-Fi, electrical disturbing devices, Microwaves; following procedure is recommended for remedies imprinting:
- Switch on the UBS 315 using key 1
 - Select laser or LED beam
 - Add in the containers the master remedies always filled in glass vials
 - Add one glass (transparent) or bottle (glass, not plastic) on the ULD 117; switch on the containers plate (key 5) and the ULD 117 for some seconds; all information are transmitted from the master substances in mini-containers into the liquid located on ULD 117. This operation takes some seconds (max. 30-50 sec)
 - Like this you can improve the taste of the fluids, duplicate remedies
 - During the information swing there is a white LED blinking as control function

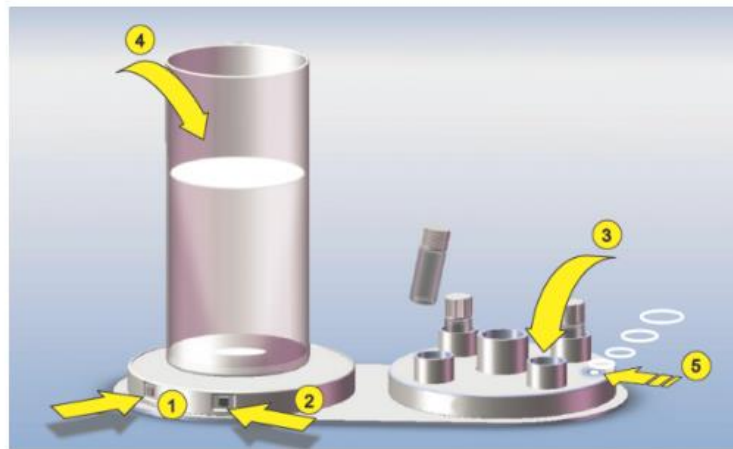

HEALTH BENEFITS:

- Cellular regeneration and deep detoxification
- Genetic modulation and genetic correction
- Pain relief, connective tissue, cartilages restoration and energizing of muscular skeletal body
- Therapy of all multiple chronic conditions
- Improvement of rheological properties of body fluids, blood, lymph
- Mitochondrial activation, energy boost
- Stress relief, depression, anxiety and insomnia therapy, emotional trauma cure
- Hormonal, metabolic and self-healing stimulation
- Cellular regeneration especially after injuries, surgeries, cell loss
- Anti-inflammatory
- Brain stimulation and hemispheres synchronization

HOW TO APPLY THE THERAPY:

- **Nervous system:** Solar plexus, Back of the head, Vertebra C7
- **Immune system:** Thymus, On and below the navel, Back of the head
- **Endocrine system:** Top of the skull, thyroid, pubic region, center of the spine
- **Respiratory system:** on the middle of sternum, under the collar bones, left and right joints shoulder thorax, on the back side between shoulder blades
- **Circulatory system:** on each kidney, underneath vertebra C7, back of the head
- **Detox and Lymphatic system:** on vertebra C7, under collar bones, back of the head, liver, each kidney
- **Gastro-intestinal tract:** Solar plexus, on the navel and below navel 1 inch, left and right of the knees, back of the head, underneath shoulder blades left and right
- **Kidneys and Urinary Tract:** back of the head, each kidney, inside the legs on the ankles, bladder and along the spine left and right (bladder meridian)
- **Liver and Gallbladder:** underneath right ribcage, back of the head, right shoulder
- **Joints and spine:** treat always the left and the right knees, shoulders, elbows, feet, hand wrists; treat inside and outside of the joints, on the top and underneath the joints; for the spine insist on vertebra C7 and sacrum, go along the spine up and down, always left and right of the spine not on the spine directly
- **Jaws, teeth, mouth and gum:** on the sinuses forehead and left and right the nose, underneath jaws, back head, behind ears and on temporal lobes
- **Eyes:** **use only the LED light, never use the laser radiation:** close the eyes and apply on each eye, underneath eyes, back of the head, temporal lobes, forehead
- **Pain:** apply directly on the painful zone(s) or on pain trigger points indicated by your health care giver
- **Wounds/Injuries:** Never apply directly on the wounds, burns, only around them
- **Skin eczemas:** Never apply directly on the wounds, burns, only around them
- **You can apply ULD 117 on acupuncture points, reflex zones**

MAINTENANCE AND SAFETY:

1. Both the laser diodes and LED have a lifespan of 20000 hours, corresponding to a usage of 4 hours per day, it means the device is lasting 10 years
2. The device enclosure is biological compatible; never place the device in humid spaces or under the direct sunlight, do not put objects heavier than 0.5 grams
3. Avoid that the small laser windows are getting dirty; clean the enclosure with a soft cotton cloth slightly damp, do not use solvents; pay attention that no water is penetrating inside of the enclosure
4. Always charge the battery, never store the device with empty battery; charging time is 2-4 hours